**HuT Labs**  
Humanitarian Technology Labs

# NETRAVAAD

— TRAINING MANUAL

# **TABLE OF CONTENTS**

- Introduction
- Package contents
- Device Overview
- Powering ON the Device
- Getting Started
- User Interface
- Basic Operation
- Maintenance
- Troubleshooting
- Safety and Precautions
- Technical Specifications
- Warranty and Support

# **1. INTRODUCTION**

The main goal of Netravaad is to enable people with speech impairment to communicate with the world using simple eye signs. The software aims to provide accurate and efficient eye sign recognition for seamless communication. This user manual provides step-by-step instructions on how to set up and use the device effectively.

## **2. PACKAGE CONTENTS**

- a. Netravaad Device
- b. Power cable
- c. Camera
- d. HDMI Cable
- e. User Manual
- f. Mini PC

## **3. DEVICE OVERVIEW**

Netravaad consists of the following components:

- a. Front-facing camera
- b. Touch Screen
- c. On/Off switch

## **4. POWERING ON THE DEVICE**

To turn on the device:

Connect the cable and switch on the device.

After powering on, the GUI will upstart, displaying the home page with four buttons: Start, User Manual, Settings, and Power Off.

## **5. GETTING STARTED**

Powering ON/OFF:


- To turn the device on, use the flip button on the side of the box. Flip it up to turn the device on and flip it down to turn the device off.
- To power OFF the device, select "Power OFF" button shown in the GUI Home screen. After touching "Power OFF" button seen in the GUI home screen, a popup showing "Do you want to Turn OFF device?" will appear on the screen. If the user clicks on the option "Yes" the device will turn OFF and if clicks "No" the device will not turn OFF.

## **6. USER INTERFACE**

The User Interface (UI) in Netravaad plays a critical role in enabling users to interact with the system effectively and efficiently. The UI is designed to facilitate communication between the user and the mobile application, making it easy for users to input their intentions and receive meaningful outputs. Given the nature of Netravaad as an assistive communication system based on eye movement recognition, the UI needs to be intuitive and accessible to users with diverse abilities.

#### A. Home Page:

The Home page of the GUI consists of four selection buttons such as “Start”, “User Guide”, “Settings” and “Power off” buttons and a Wi-Fi connection indication screen as shown in the below Figure 1.

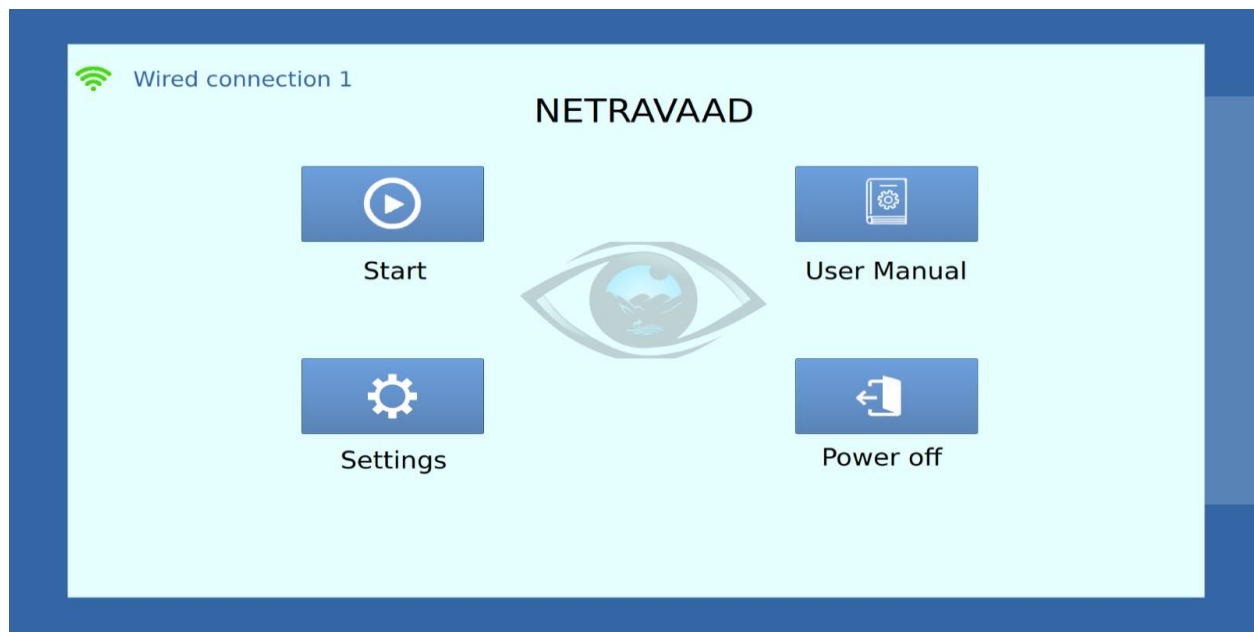**Fig. 1 Home screen page**

#### B. Settings Page, User Manual Page, Power off


a. Settings page consist of voice prompt, volume control, Instruction language selection and Instruction display turn off as shown in Fig. 2.

b. Voice Prompt and Volume Control:

- Turn ON/OFF voice prompt to hear instructions.
- Adjust volume using the volume control button.


c. Language Mode and Instruction Display:

- Select between English, Malayalam, Hindi for instructions.
- Turn ON/OFF the Instruction Display button.

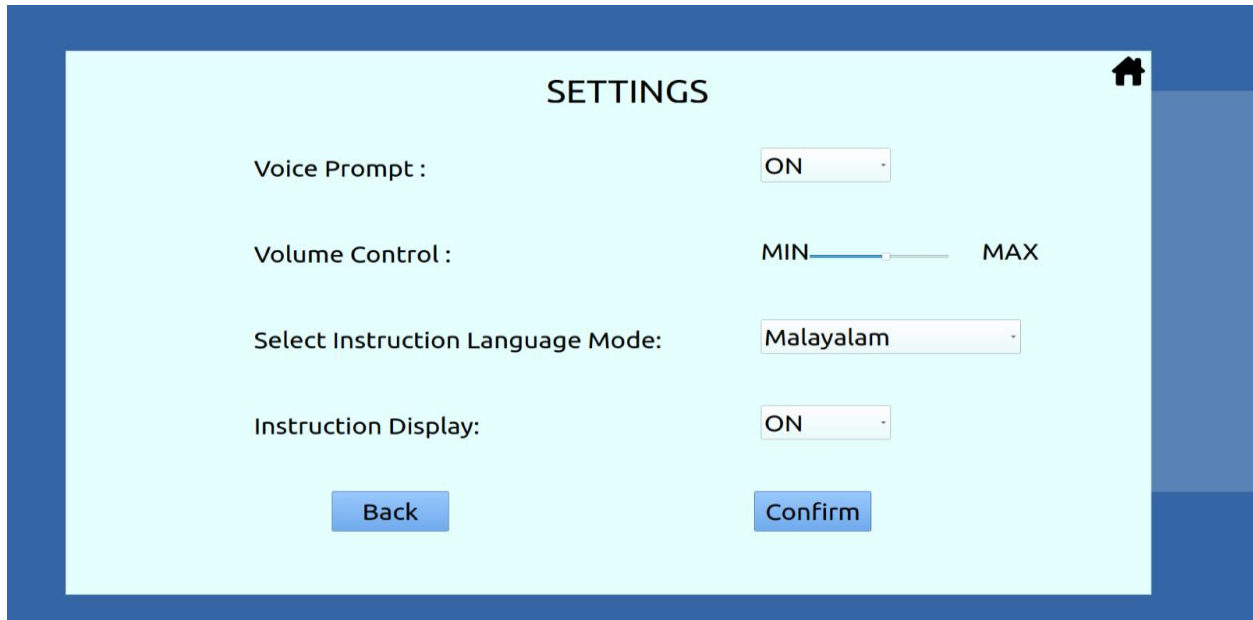**Fig. 2 Settings Page**


d. User Manual button will direct you to the user manual page in Fig.3

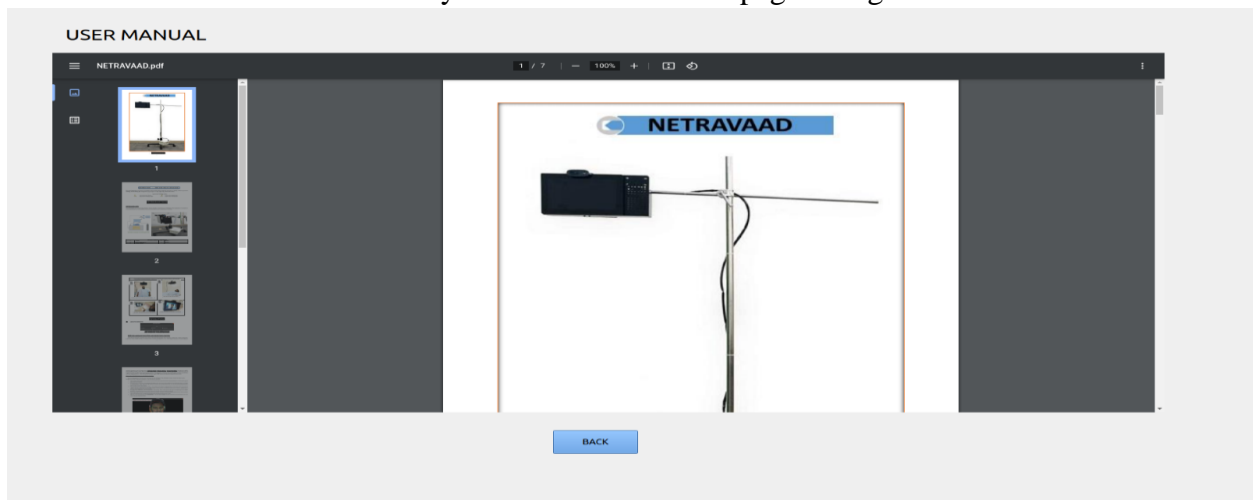**Fig. 3 User Manual**

e. Power off button- After touching “Power OFF” button seen in the GUI home screen, a popup showing “Do you want to Turn OFF device?” will appear on the screen. If the user clicks on the option “Yes” the device will turn OFF and if clicks “No” the device will not turn OFF.


Clicking the start button will be directing you to the calibration page as shown in Fig.4

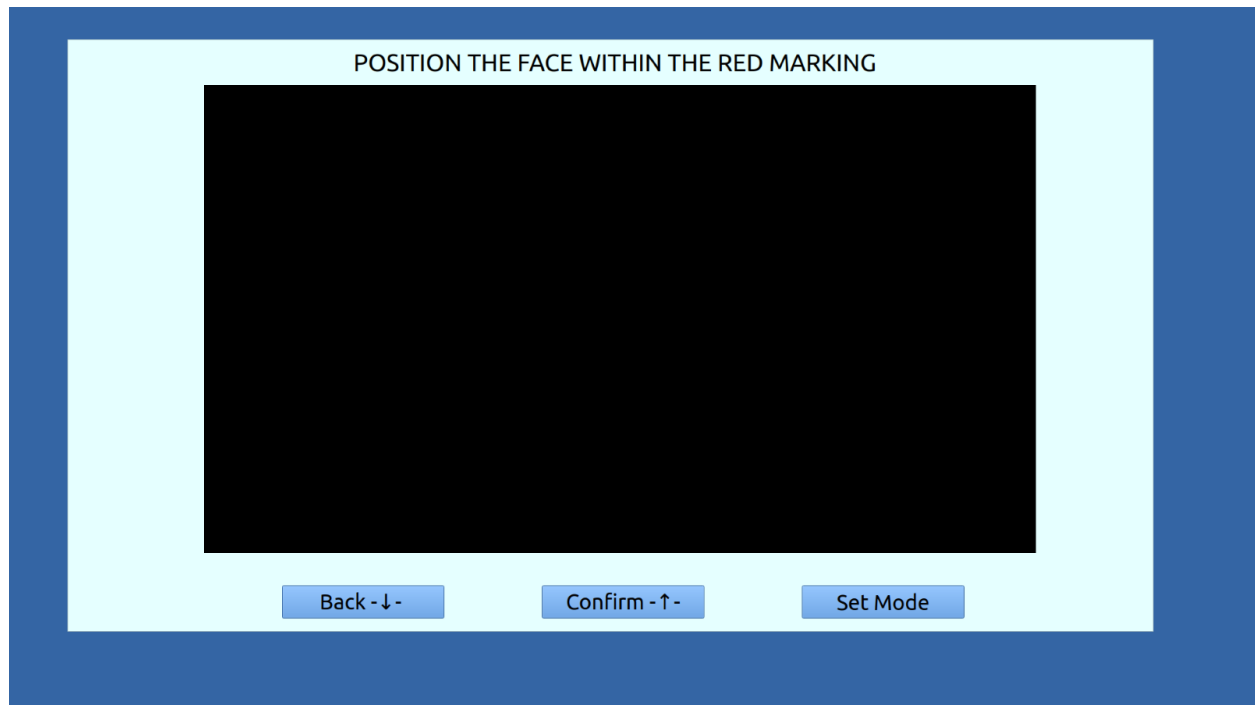**Fig.4 Calibration Page**

### **C. Start Button:**


- **Head Position Correction and Instruction charts:**
  - Click the "Start" button to navigate to the head position correction page.
  - Correct your head position based on the displayed image.
  - Follow the instruction charts and click the skip button if you already know the instructions and eye charts.
- **Confirming Head Position:**
  - Click the "Confirm" button or use eye patterns to confirm.
- **Back and Set Mode Buttons:**
  - Use the "Back" button to return to the home page.
  - Click "Set Mode" to go to the mode selection page as shown in Figure 5 for mode selection.

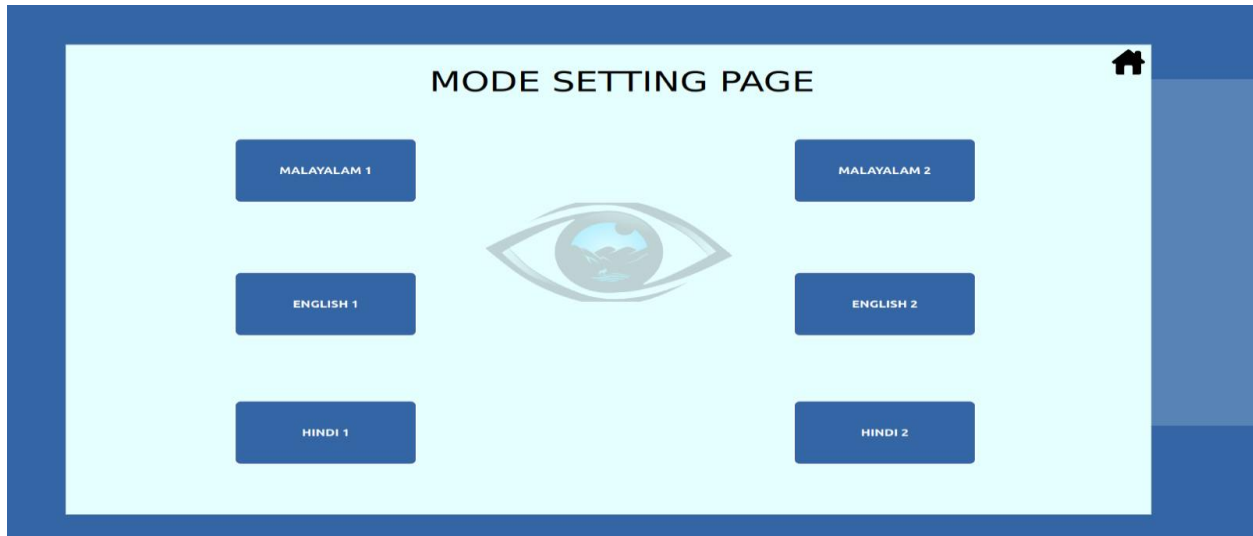**Fig. 5 Mode Setting Page**

#### **D. Mode Selection:**


Users can choose the desired mode of operation (English 1, English 2, Malayalam1, malayalam2, Hindi1, Hindi 2) from the main screen. This could be achieved by showing eye signs from the user. After clicking the confirm button or using eye patterns for confirm, you will be directed to the mode selection pages. If the user looks Straight-Close-Straight eye pattern, then it will move to the next mode and to enter inside the mode the user need to look Straight-Top-Straight. The below Figure 6 shows the sample of mode selection page.

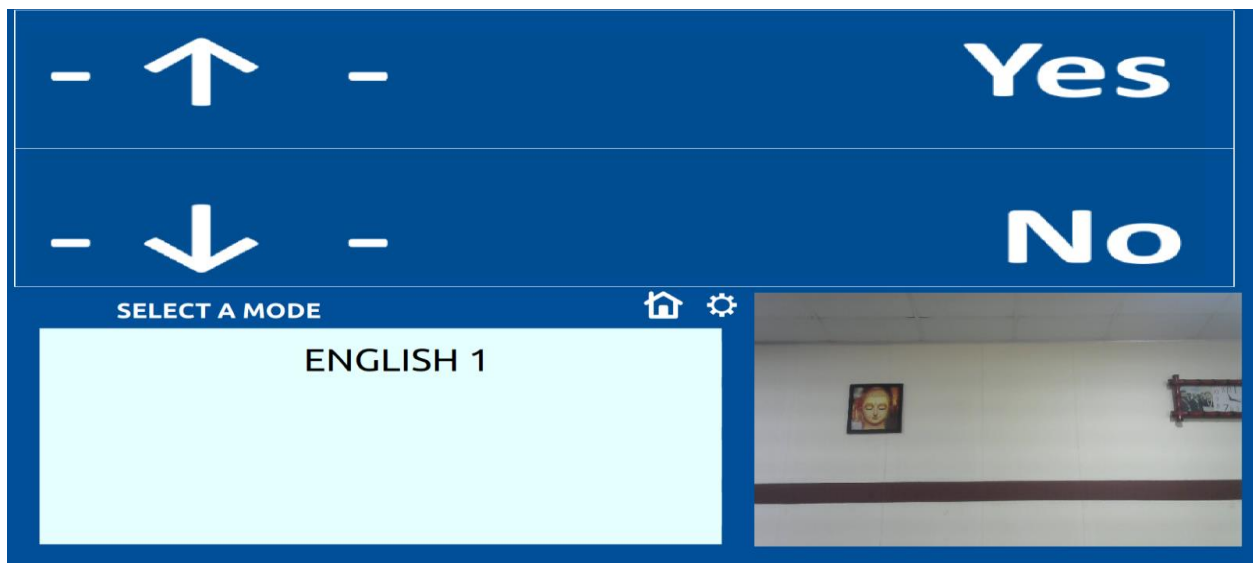**Fig. 6 Mode Selection Page**

After clicking the confirm button, you will be generated to the English 1 mode.

**English 1** mode consists of predefined words shown in Fig.7 You can select the suitable word you wanted

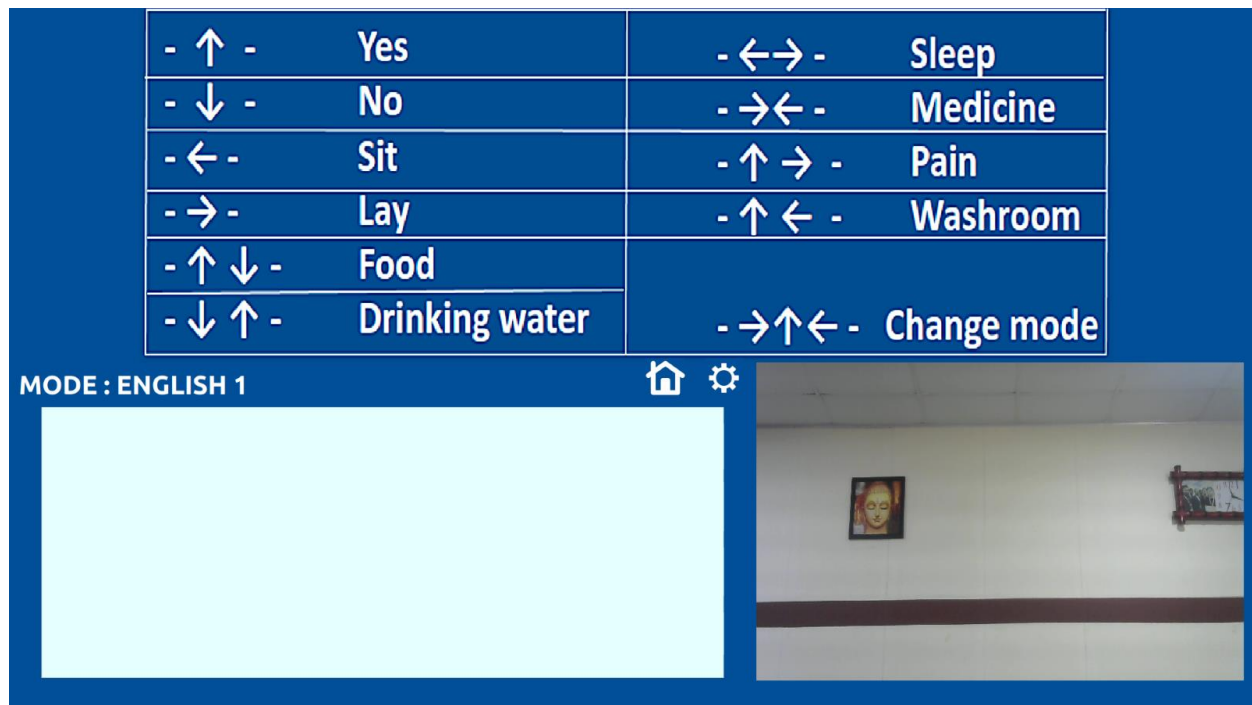**Fig. 7 English Mode 1**

For example- If you want to select the word “Sit” -Use the eye pattern-

First look straight with your eyes, then look left with your eyes, again look straight with your eyes

The selected word will be displayed in the screen and a voice output for the selected word will be generated

Similarly, you can select the suitable word you wanted from the chart

10-For moving to the next mode English 2, Use the change mode pattern

“Change Mode” - Look straight with your eyes

Look right with your eyes


Look Up with your eyes

Look left with your eyes

Look straight with your eyes

If you had done the correct eye pattern, you will be generated to the English 2 Mode shown in Figure 8.

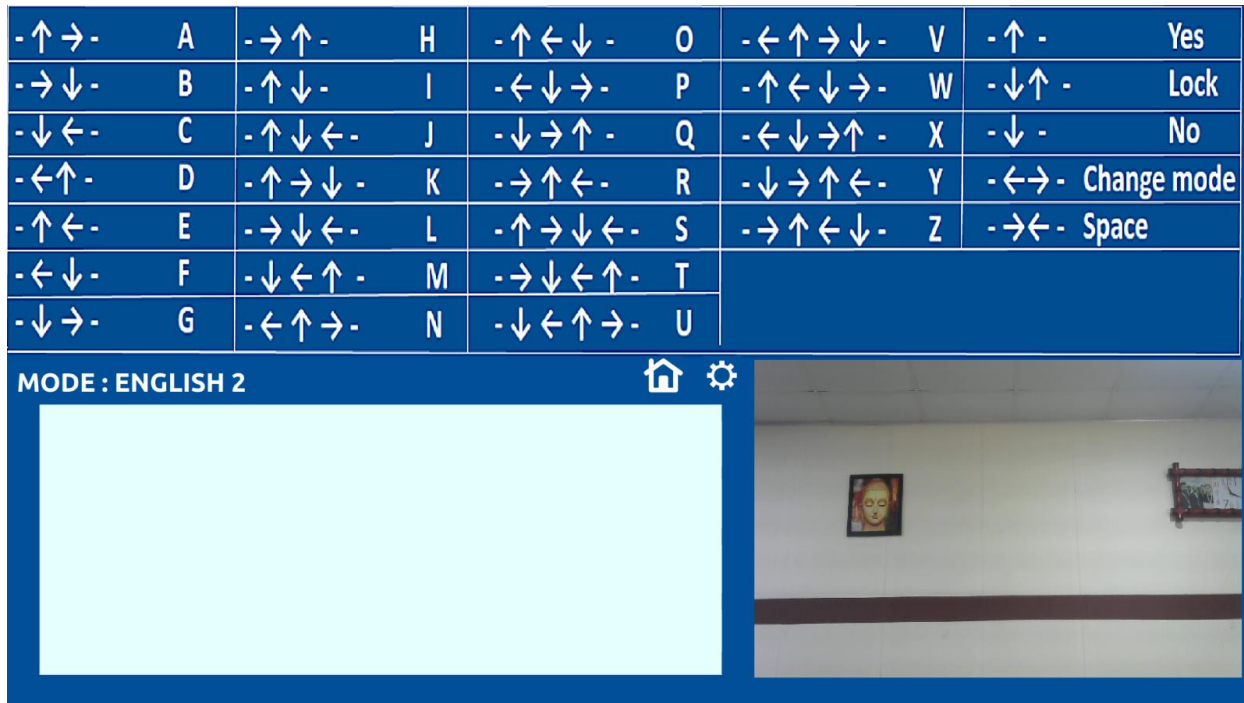**Fig. 8 English Mode 2**

**English Mode 2** consist of 26 alphabets in English with 26 eye patterns

First for word prediction, select an alphabet

For example, select the letter “a” by using the eye pattern

Look straight with your eyes

Look Up with your eyes

Look Right with your eyes

Look straight with your eyes

By selecting the letter “a “, it will be displayed in the screen and use the “yes” pattern to confirm the letter

After confirming the letter, a word starting with the letter “a” will be displayed, use the eye pattern “no” for moving to the next word starting with letter “a”

Likewise, till it shows the suggestions over, you can see the different words starting with letter “a “

If you want to select a word, then use the eye pattern “yes”

Similarly, you can select any letter you want and the corresponding words starting with that letter.

For generating numbers, select the letter “N”


It will display the text “Words starting with N”

Use the eye pattern for “NO”

It will display the text “Numeric word”

Use the eye pattern for “YES”, it will generate you to numeric mode

Numeric mode consists of numbers from 0 to 9 as shown in figure 9

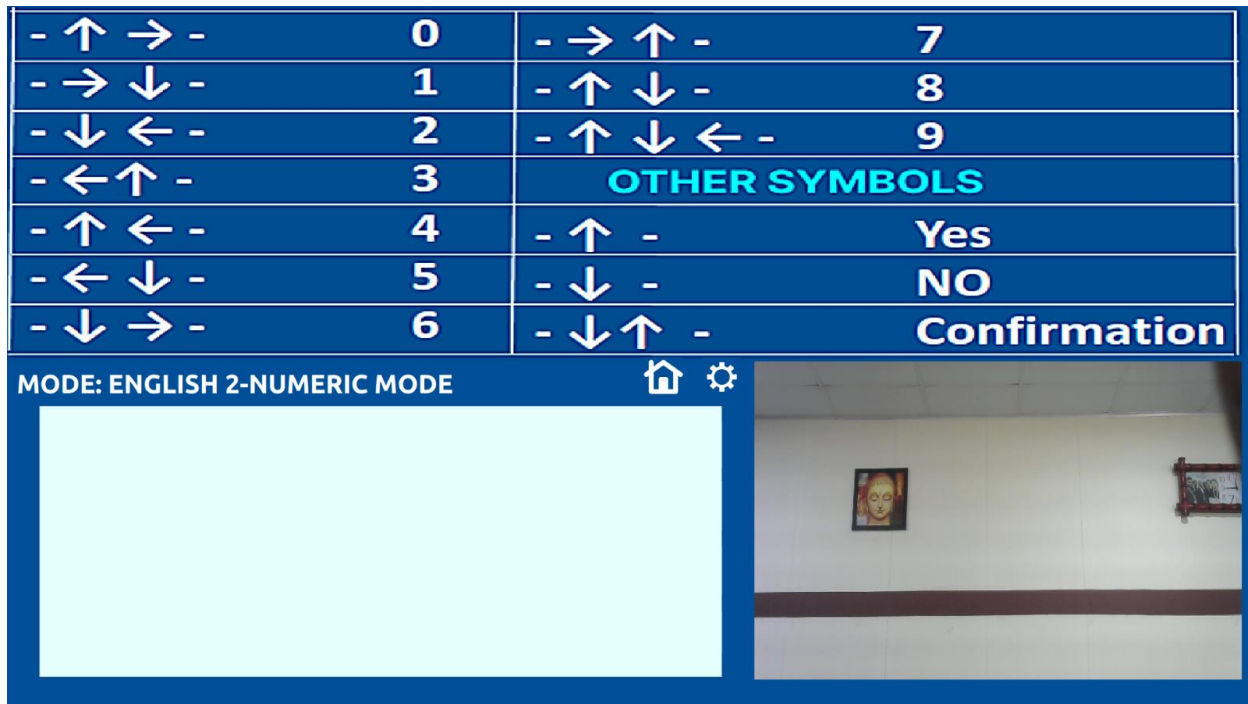**Fig. 9 Numeric Mode**

For selecting a number example “2” Use the eye pattern provided for number “2”

Then use the eye pattern provided for confirmation to confirm the number selected.


For generating the sentences, Select the letter “S”,

It will display the text “Words starting with S”

Use the eye pattern for “NO”

It will display the text “Sentence Mode”

Use the eye pattern for “YES”, it will generate you to sentence mode as shown in figure 10

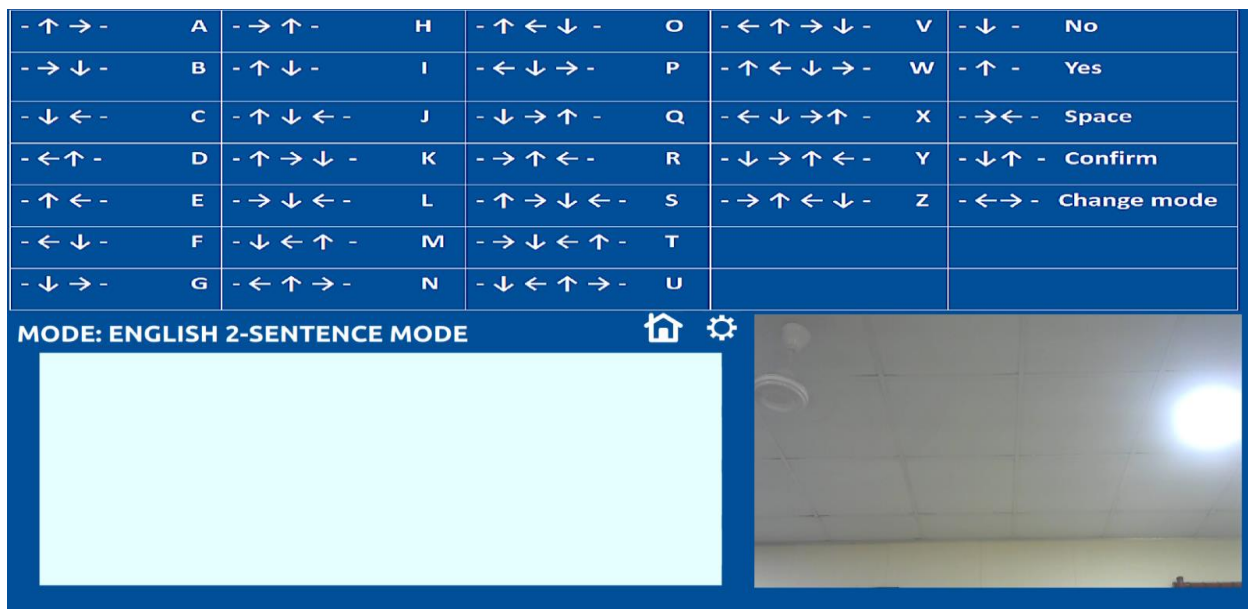
**Fig. 10 Sentence Mode**

Here you can generate sentences you wanted

Example - “I need food”


First select the letter “I” then use the space eye pattern for providing space after the letter “I”, then using the eye patterns select the letters “N”, “E”, “E”, “D”, then use the eye pattern for space and at last select the letters “F”, “O”, “O”, “D” using the eye patterns

For each selecting letter, 3-word suggestions will be shown as shown in figure 11, if you want any of these words, select that or ignore.

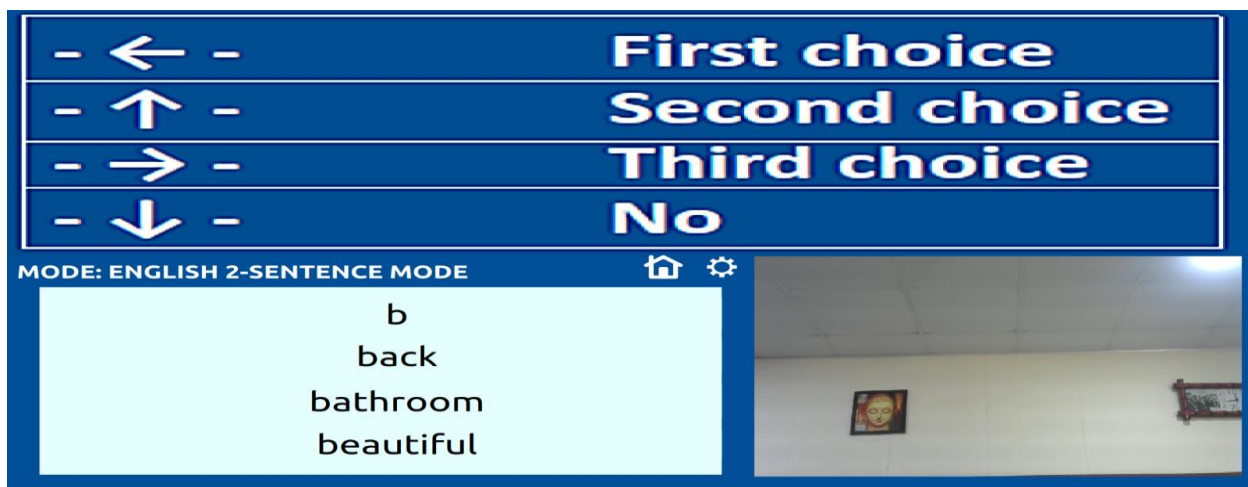**Fig. 11 Word suggestion**

To confirm the sentence, use the eye pattern to confirm (look straight with your eyes, close your eyes, look up with your eyes, look straight with your eyes), it will generate you to the confirmation page as shown in Fig.12.

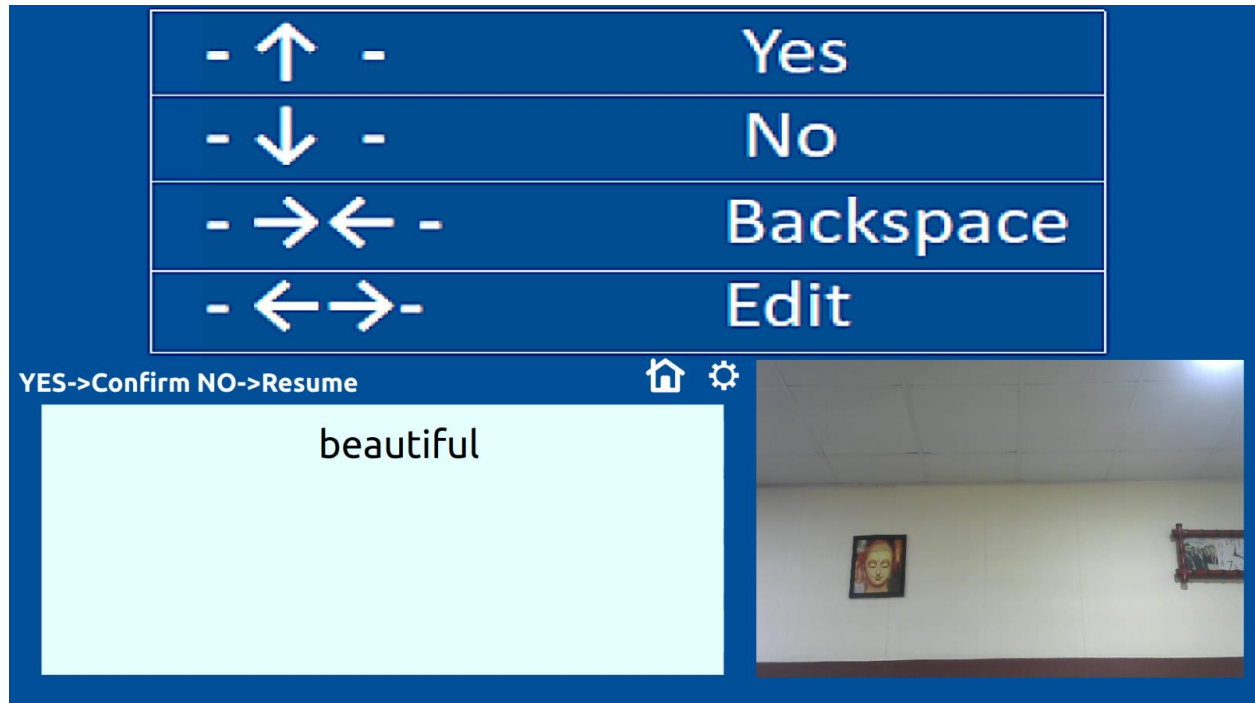**Fig. 12 Confirmation Page**

To edit the sentences, use the eye pattern (Look straight with your eyes, look left with your eyes, look right with your eyes, look straight with your eyes) to select the option edit.

After selecting the edit option, you will be generated to the next page which consist of four options move left, move right, replace and back


To move left to the word, use the eye pattern to select the option move left and select the replace option using eye pattern to replace the letter you wanted.

To remove the letter, use the backspace option which will remove the letter till you wanted

Use the confirm option to confirm the modified sentence

Next, to move to the next mode, use the change mode pattern, which will generate you to the **Malayalam Mode 1**

Malayalam 1 mode consists of predefined Malayalam words, you can select the suitable word you wanted as shown in figure 13.

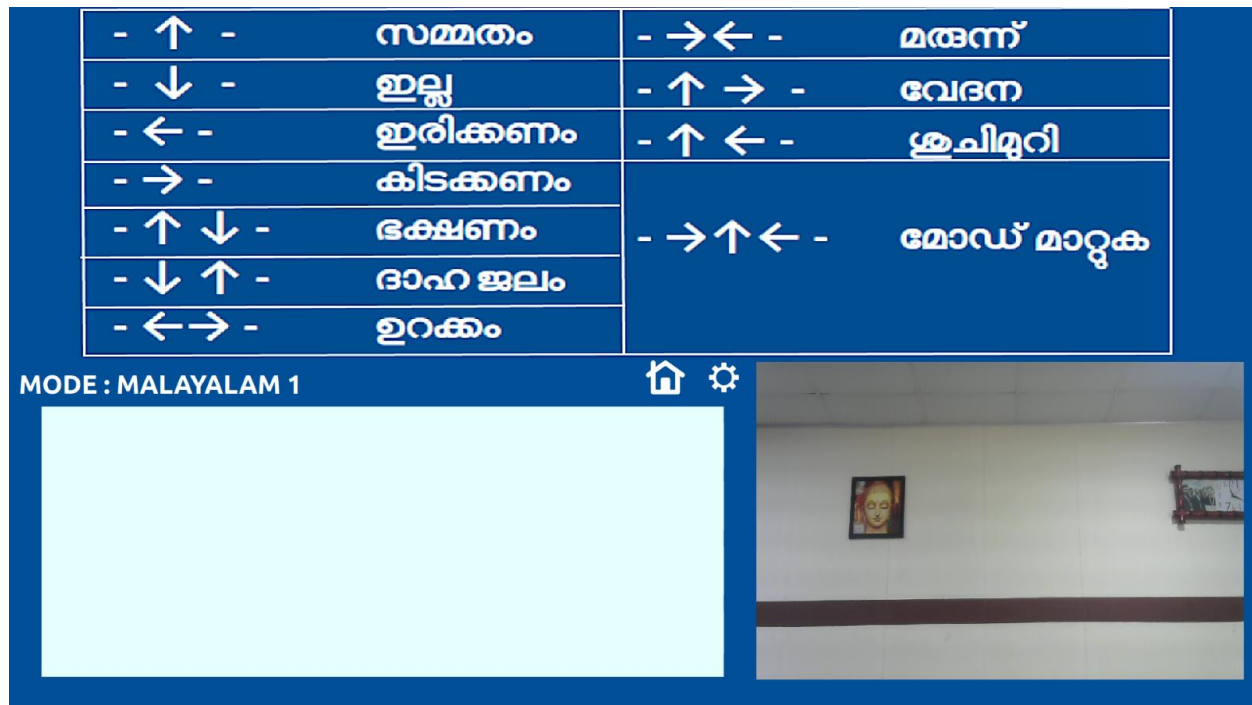**Fig. 13 Malayalam Mode 1**

For example, if you want to select the word “ഇരിക്കണം” use the eye pattern:

First look straight with your eyes, then look left with your eyes, again look straight with your eyes

The selected word will be displayed in the screen and a voice output for the selected word will be generated

Similarly, you can select the suitable word you wanted from the chart

12-For moving to the next mode Malayalam 2, Use the change mode pattern

“Change Mode” - Look straight with your eyes

Look right with your eyes


Look Up with your eyes

Look left with your eyes

Look straight with your eyes

If you had done the correct eye pattern, you will be generated to the Malayalam 2 Mode

Malayalam Mode 2 consist of 29 selected letters in Malayalam with 29 eye patterns as shown in Figure 14

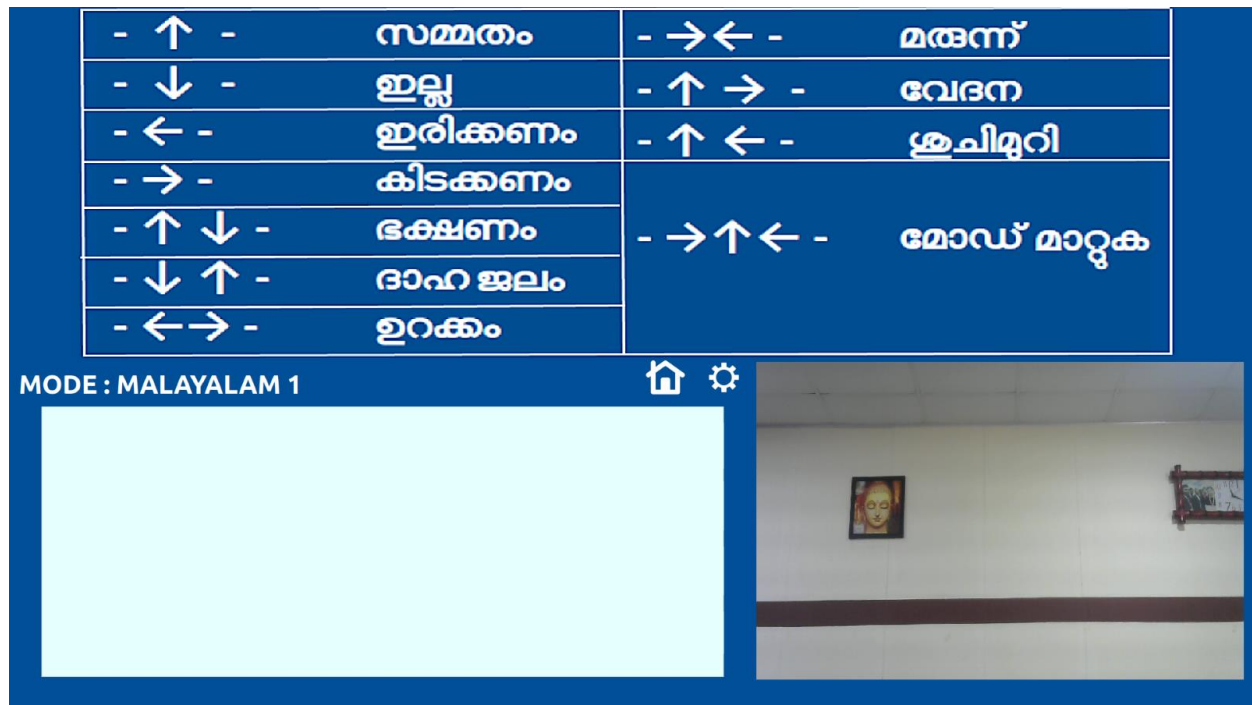**Fig. 13 Malayalam Mode 1**

For example- If you want to select the word “ഇരിക്കണം” -Use the eye pattern-

First look straight with your eyes, then look left with your eyes, again look straight with your eyes

The selected word will be displayed in the screen and a voice output for the selected word will be generated

Similarly, you can select the suitable word you wanted from the chart

12-For moving to the next mode Malayalam 2, Use the change mode pattern

“Change Mode” - Look straight with your eyes

Look right with your eyes

Look Up with your eyes


Look left with your eyes

Look straight with your eyes

If you had done the correct eye pattern, you will be generated to the Malayalam 2 Mode

**Malayalam Mode 2** consists of 29 selected letters in Malayalam with 29 eye patterns as shown in figure 14.

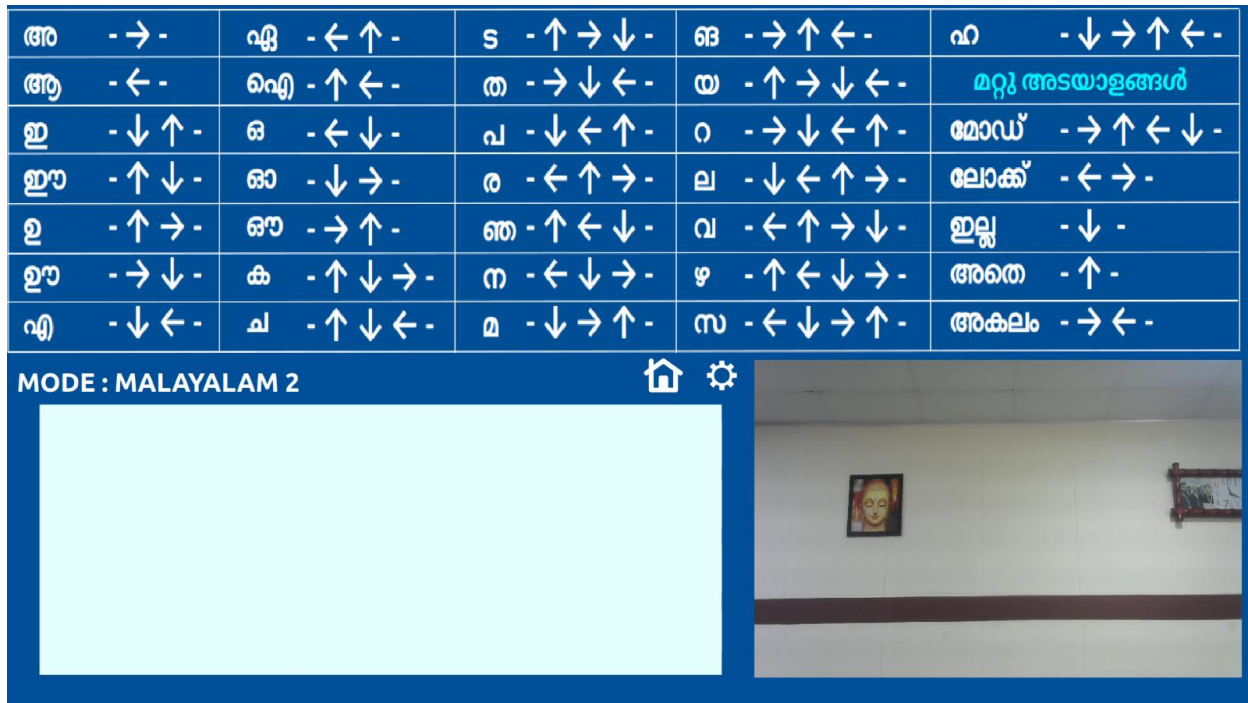**Fig. 14 Malayalam Mode 2**

First for word prediction, select one letter

For example, select the letter “അ” by using the eye pattern

Look straight with your eyes

Look Right with your eyes

Look straight with your eyes

By selecting the letter “അ”, it will be displayed in the screen and use the “അതെ” pattern to confirm the letter

After confirming the letter, a word starting with the letter “അ” will be displayed, use the eye pattern “ഇല്ല” for moving to the next word starting with letter “അ”

Likewise, till it shows the “നിർദ്ദേശം കഴിഞ്ഞു”, you can see the different words starting with letter “അ”

If you want to select a word, then use the eye pattern “അതെ”

Similarly, you can select any letter you want and the corresponding words starting with that letter.

For generating numbers, select the letter “സ”

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
Use the eye pattern for “അതെ”, it will generate you to numeric mode

Numeric mode consists of numbers from 0 to 9 as shown in figure 15

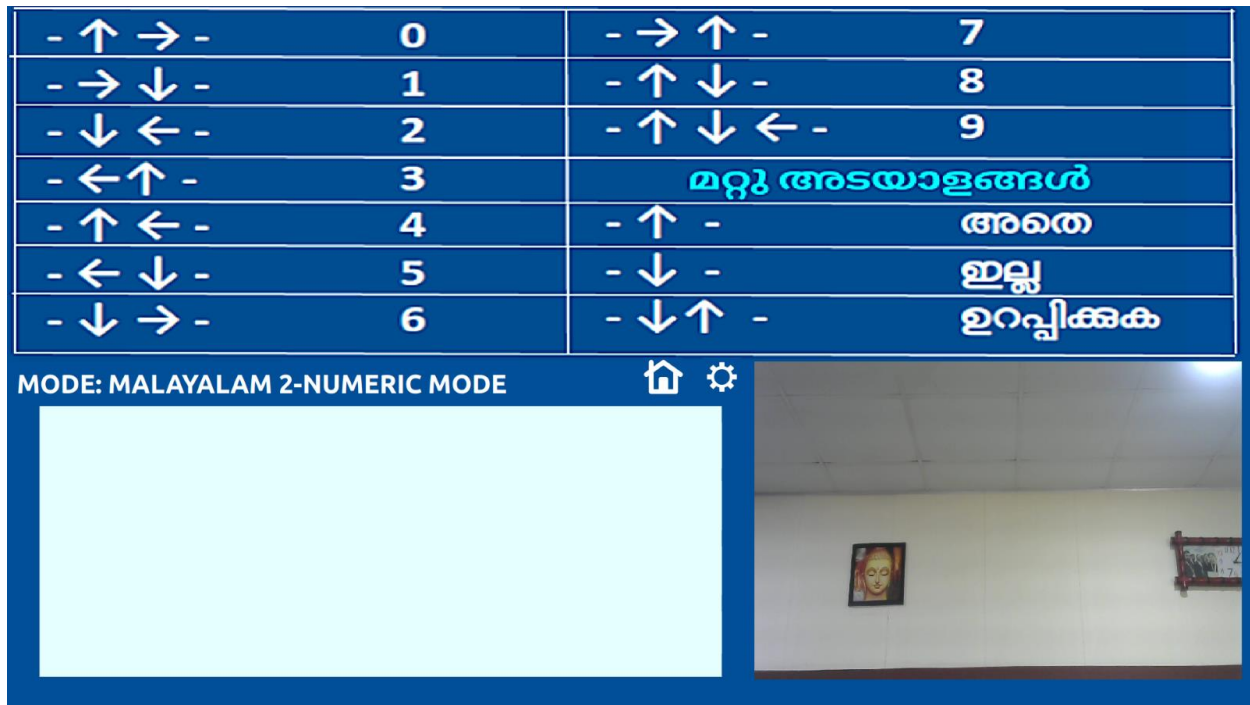**Fig. 15 Malayalam Numeric Mode**

For selecting a number example “2” Use the eye pattern provided for number “2”

Then use the eye pattern provided for confirmation to confirm the number selected.

The selected number with Malayalam audio output will be generated.

For generating the sentences, Select the letter “വ”,

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Use the eye pattern for “അതെ”, it will generate you to sentence mode as shown in figure 16.

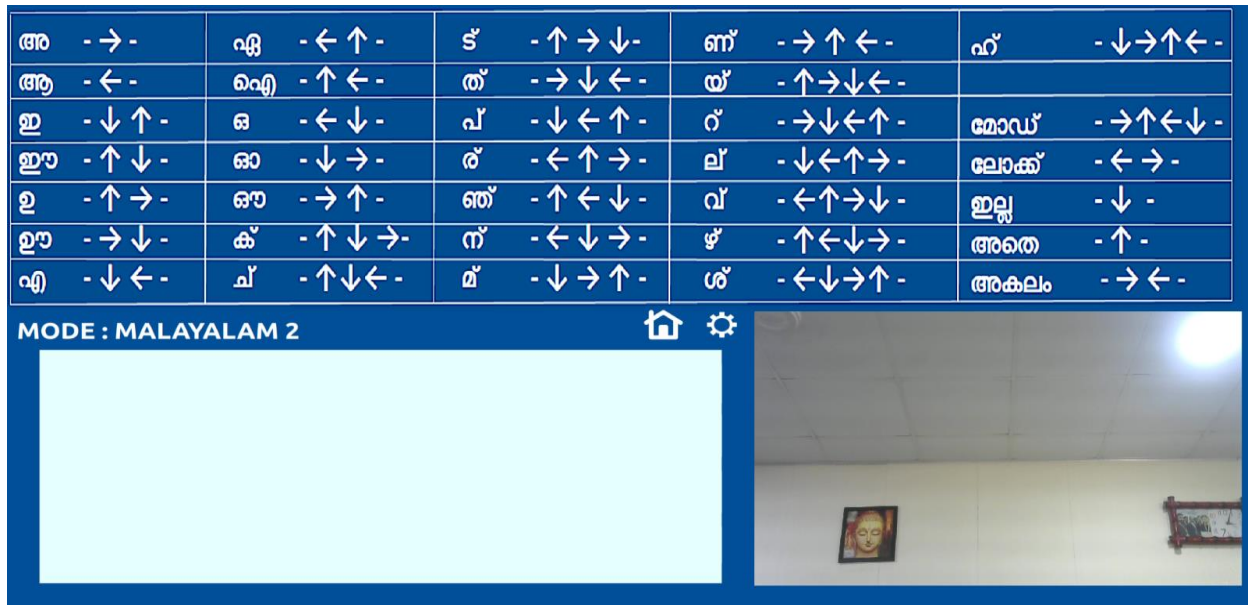
Fig. 16 Sentence Mode

Here you can generate sentences you wanted

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First select the letter “എ” then select the letter “ന്”, then select the letter “ഇ”, then select the letter “ക്”, then again select the letter “ക്”, you will now generate the word “എനിക്കു്”.

For putting space use the eye pattern for “അകലം”

For each selecting letter 3-word suggestions will be shown as shown in figure 17, if you want any of these words, select that or ignore.


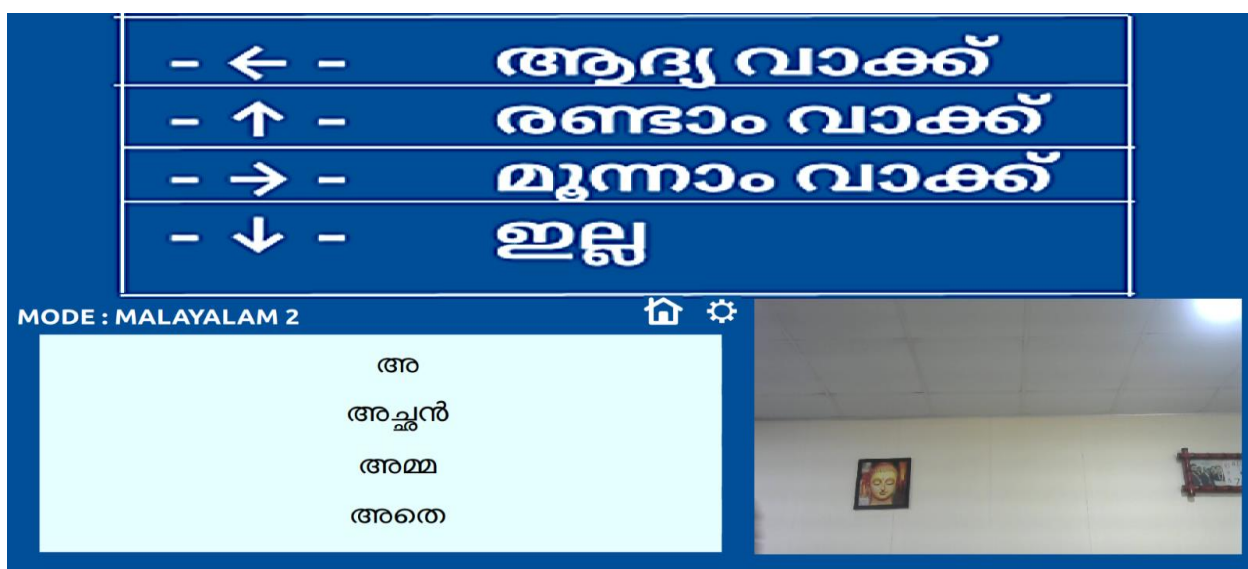
Fig. 17 Malayalam Word suggestion

To confirm the sentence, use the eye pattern to "ലോക്ക്" (Look straight with your eyes, look left with your eyes, look right with your eyes, look straight with your eyes), it will generate you to the confirmation page as shown in figure 18.

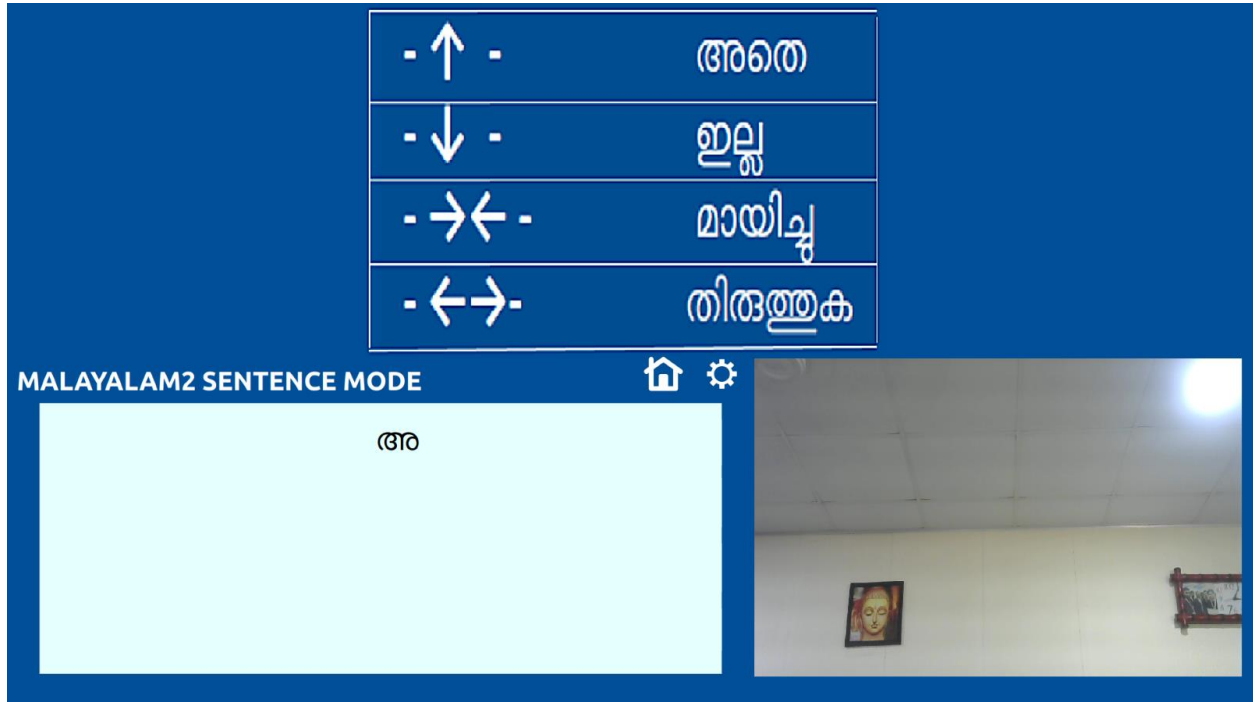**Fig. 18 Malayalam Lock Page**

To edit the sentences, use the eye pattern for “തിരുത്തുക” (Look straight with your eyes, look left with your eyes, look right with your eyes, look straight with your eyes) to select the option “തിരുത്തുക.”

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
To remove the letter, use the “മായിച്ചു” option which will remove the letter till you wanted

Use the "ലോക്ക്" option to confirm the modified sentence

To use screen on/off mode, use the change mode option in Malayalam mode 2, which will generate you to the Hindi Mode1, Hindi Mode2 then screen on/off mode selecting that will make the screen blank and to resume back use eye pattern displayed in the screen blank screen.

Next, to move to the next mode, use the change mode pattern, which will generate you to the **Hindi Mode 1**

Hindi 1 mode consist of predefined Hindi words as shown in figure 19, you can select the suitable word you wanted

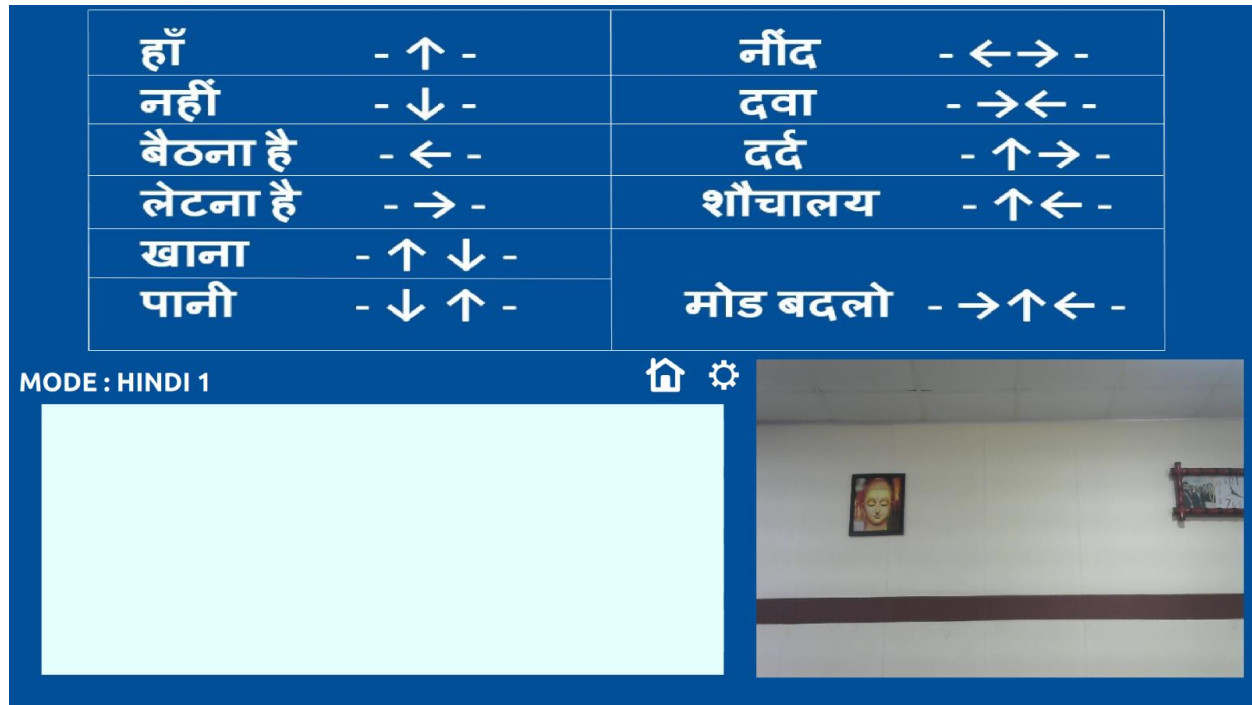**Fig. 19 Hindi Mode 1**

For example- If you want to select the word “बैठना है” - Use the eye pattern-

First look straight with your eyes, then look left with your eyes, again look straight with your eyes

The selected word will be displayed in the screen and a voice output for the selected word will be generated

Similarly, you can select the suitable word you wanted from the chart

For moving to the next mode Hindi Mode 2, Use the change mode pattern

“मोड बदलो” - Look straight with your eyes

Look right with your eyes

Look Up with your eyes

Look left with your eyes

Look straight with your eyes


If you had done the correct eye pattern, you will be generated to the Hindi 2 Mode


**Hindi Mode 2** consists of 29 selected letters in Hindi with 29 eye patterns as shown in figure 20.

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MODE : HINDI 2

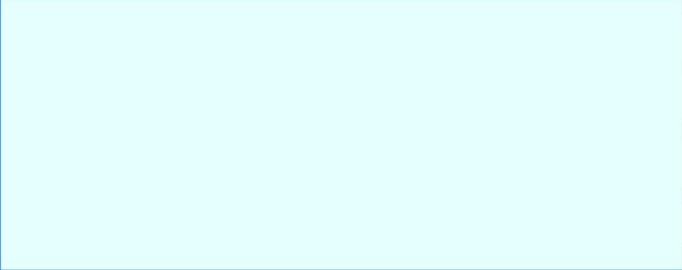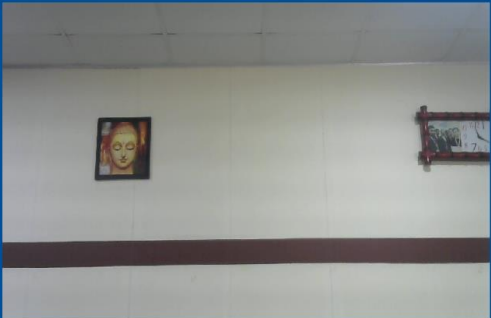**Fig. 20 Hindi Mode 2**

First for word prediction, select one letter

For example, select the letter “अ” by using the eye pattern

Look straight with your eyes

Look Right with your eyes

Look straight with your eyes

By selecting the letter “अ”, it will be displayed in the screen and use the “हाँ” pattern to confirm the letter

After confirming the letter, a word starting with the letter “अ” will be displayed, use the eye pattern “नहीं” for moving to the next word starting with letter “अ”

Likewise, till it shows the “कोई निर्देश उपलब्ध नहीं है”, you can see the different words starting with letter “अ”

If you want to select a word, then use the eye pattern “हाँ”

Similarly, you can select any letter you want and the corresponding words starting with that letter.

For generating numbers, select the letter “स”

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
It will display the text “हिंदी अंक”

Use the eye pattern for “हाँ ”, it will generate you to numeric mode

Numeric mode consists of numbers from 0 to 9 as shown in figure 21.

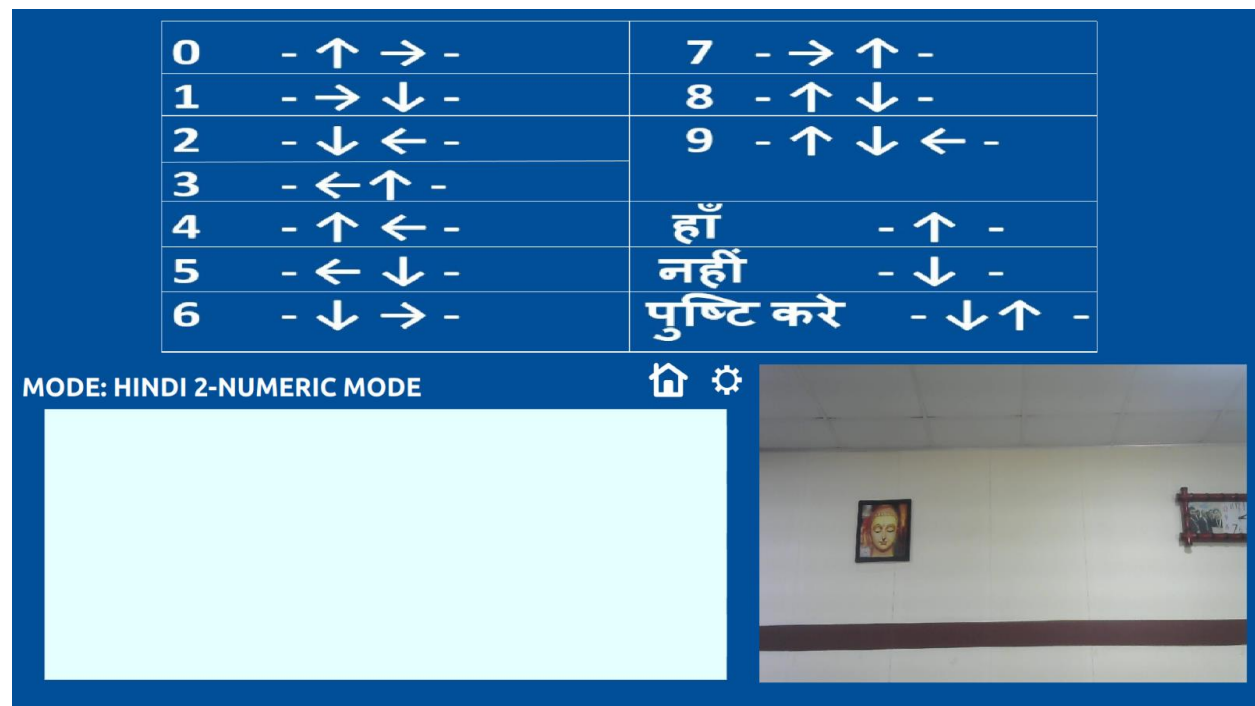**Fig. 21 Hindi Numeric Mode**

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Then use the eye pattern provided for “पुष्टि करे” to confirm the number selected.

The selected number with hindi audio output will be generated.

For generating the sentences, Select the letter “व”,

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
Use the eye pattern for “नहीं”

It will display the text “हिंदी वाचक”

Use the eye pattern for “हाँ”, it will generate you to sentence mode as shown in figure 22.

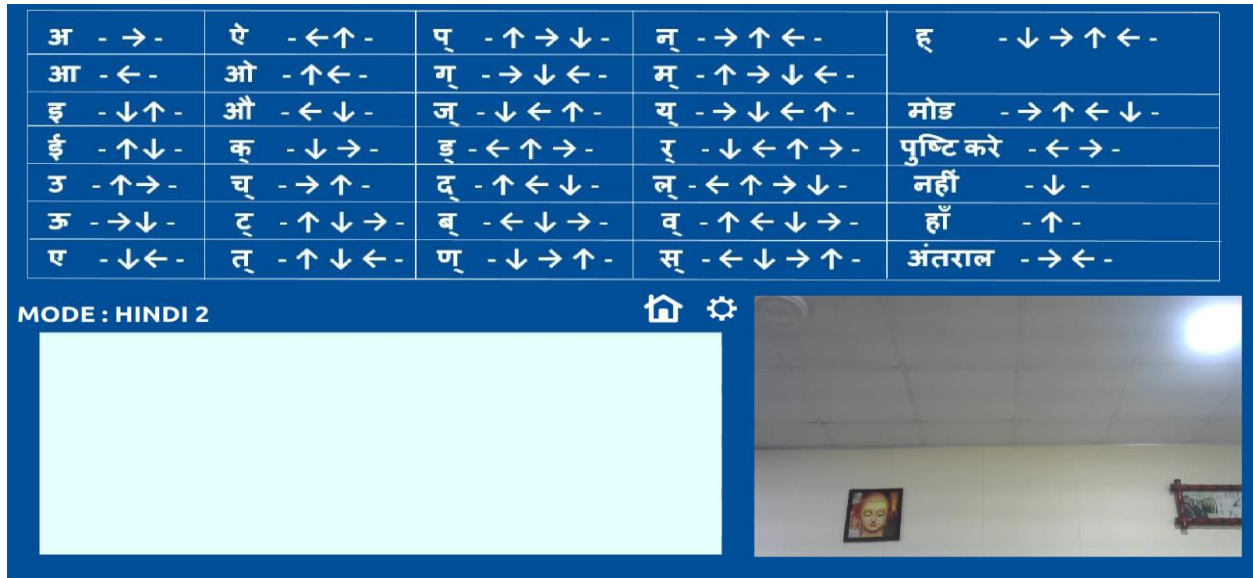**Fig. 22 Hindi Sentence Mode**

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
Example - “चाहिए”

First select the letter “च्” then select the letter "आ", then select the letter "ह्", then select the letter "इ", then again select the letter "ए", you will now generate the word “चाहिए”.

For putting space use the eye pattern for “अंतराल”

For each selecting letter 3-word suggestions will be shown in figure 23, if you want any of these words, select that or ignore by selecting “ नहीं ”.

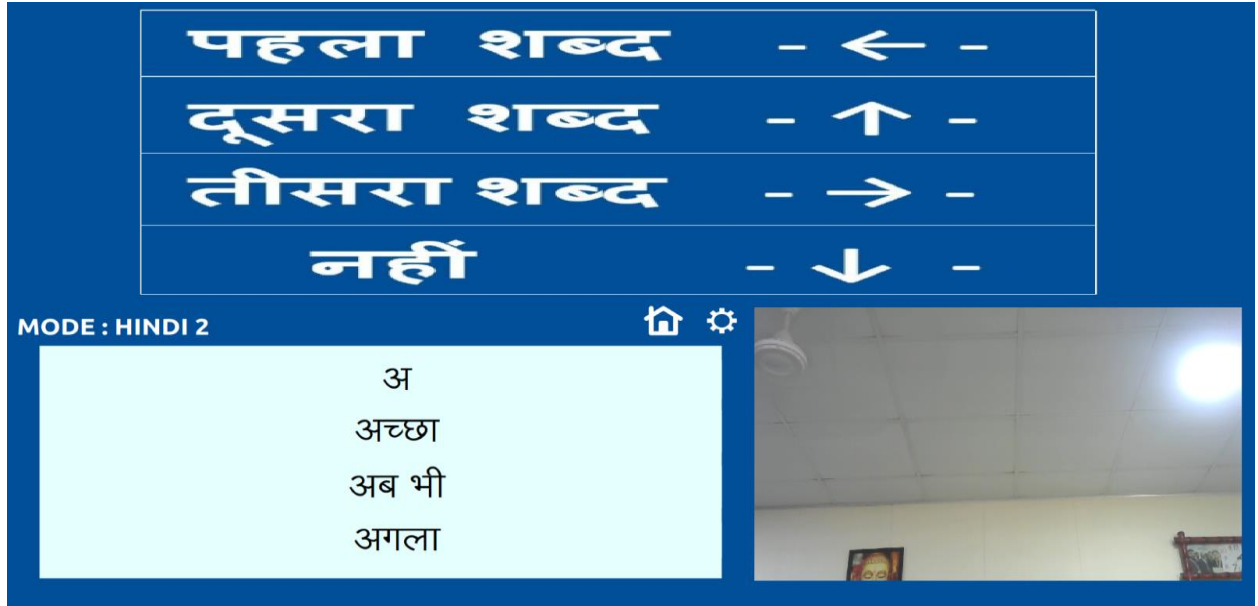
**Fig. 23 Word Suggestion**

To confirm the sentence, use the eye pattern to "पुष्टि करे" (Look straight with your eyes, look left with your eyes, look right with your eyes, look straight with your eyes), it will generate you to the confirmation page

To edit the sentences, use the eye pattern for "सुधारो" (Look straight with your eyes, look left with your eyes, look right with your eyes, look straight with your eyes) to select the option "सुधारो."

After selecting the "सुधारो." option, you will be generated to the next page which consist of four options दाईं खिसको, बाएं खिसको, बदलना, and उल्टा as shown in figure 24.

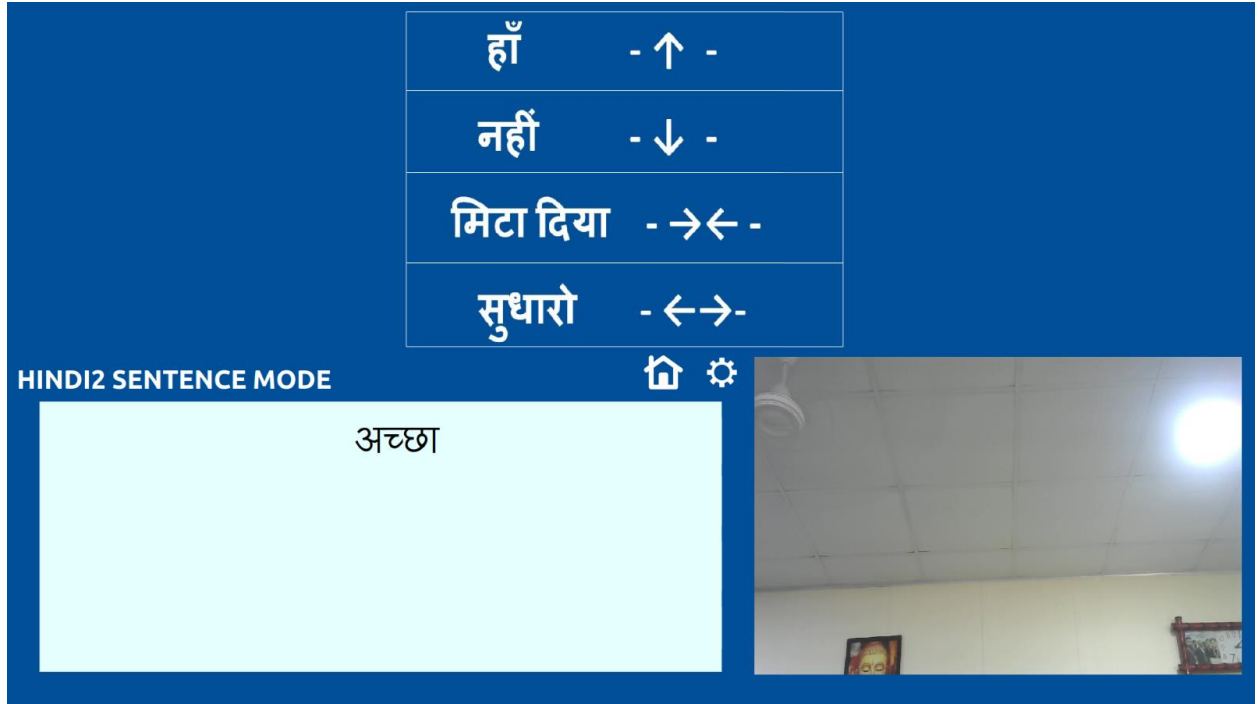
**Fig. 24 Hindi Confirmation Page**

To move left to the word, use the eye pattern to select the option “दाईं खिसको,” and select the “बदलना” option using eye pattern to replace the letter you wanted.

To remove the letter, use the “उल्टा” option which will remove the letter till you wanted

Use the “पुष्टि कर” option to confirm the modified sentence

To use screen on/off mode, use the change mode option Hindi Mode 2, which will generate you to the screen on/off mode selecting that will make the screen blank and to resume back use eye pattern displayed in the screen blank screen as shown in figure 25 and 26.

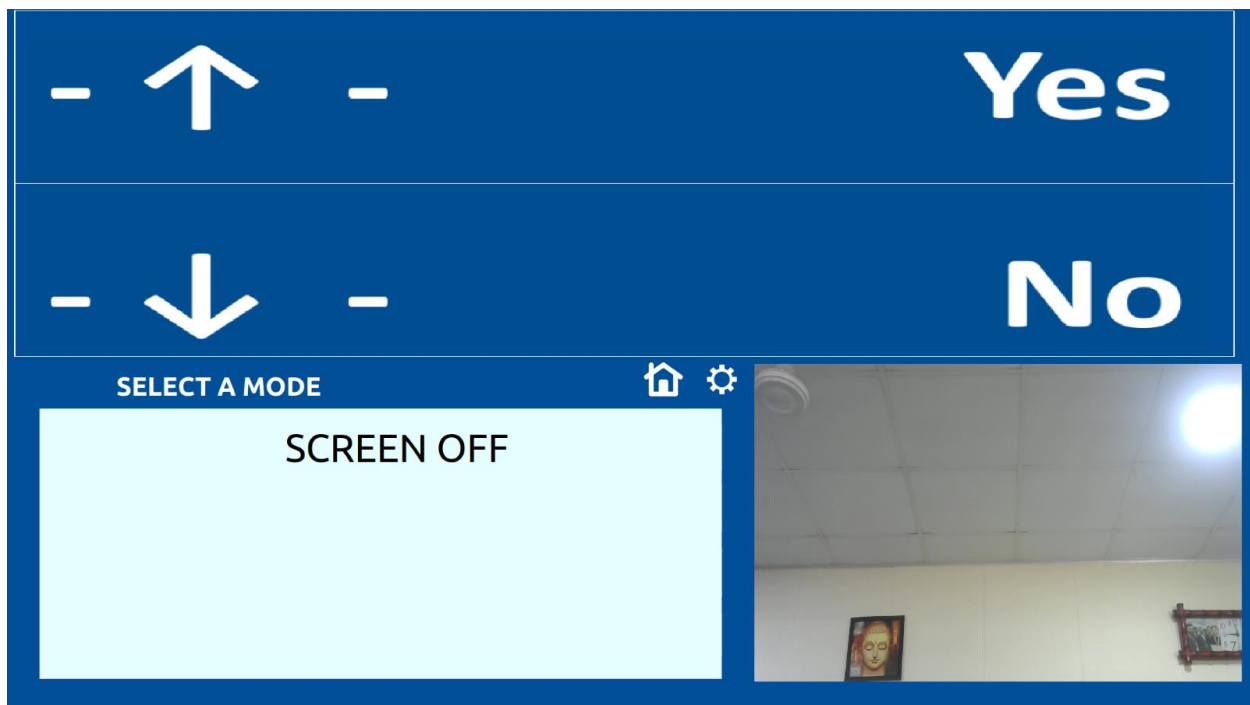**Fig. 25 Screen on/off selection**

**Fig. 26 Screen on/off page and go back pattern.**

## **7. BASIC OPERATION**

### **a. Eye Sign Detection:**

- After calibration, focus your eyes on a specific area on the screen to select an option.
- Blink or hold a gaze to confirm your choice.

### **b. Speech Output:**

- Once an option is selected, the device converts your eye signs into speech.
- Adjust the volume in the settings.

## **8. MAINTENANCE**

### **a. Cleaning:**

- Clean the device's screen and camera lens regularly with a soft, lint-free cloth.
- Avoid using harsh chemicals or abrasive materials.

## **9. TROUBLESHOOTING**

If you encounter issues, refer to the troubleshooting section in the User Manual or contact our customer support for assistance.

## **10. SAFETY AND PRECAUTIONS**

- The device is for communication purposes only.
- Do not immerse the device in water.
- Keep the device away from extreme temperatures and direct sunlight.
- Ensure it's out of reach of small children.

## **11. TECHNICAL SPECIFICATIONS**

- Screen Size: [13inch Touch Screen]
- Operating System: [UBUNTU 20.04]
- Language Support: [English & Malayalam & Hindi]

## **12. WARRANTY AND SUPPORT**

NETRAVAAD comes with a limited warranty. Please refer to the warranty card for details. For support and inquiries, contact our customer support team at [hutlabs@am.amrita.edu](mailto:hutlabs@am.amrita.edu) or <https://t2hinnovations.com/>

Thank you for choosing the Netravaad. We hope this user manual helps you make the most of this innovative communication device.

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